


Amplex UltraRed (AUR) – Apoplastic H_2O_2 Detection

2025-12-10 Ahyeon Cheon

1. Stock: 1 mg; molecular weight (MW) = 300 g/mol.
 - Stock solution: 10 mM in DMSO.
(To prepare, add approximately 333 μL DMSO directly to the original stock vial.)
2. Store at $-20\text{ }^{\circ}\text{C}$ in the dark.
3. Example use (abscission zone)
For the abscission zone shown below, AUR was diluted to 10 μM in distilled water and applied for 3 min in the dark before imaging. Because the abscission zone lacks a cuticle, AUR penetrated readily; other tissues may require different conditions.

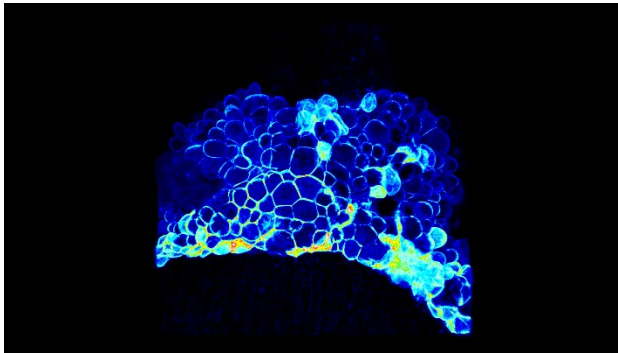

[AUR reference] Detecting hydrogen peroxide in leaves in vivo – a comparison of methods

https://onlinelibrary.wiley.com/doi/full/10.1111/j.1399-3054.2008.01176.x?utm_source=chatgpt.com